

Accord Medical Products Pvt. Ltd.



**Manufactures Of : Physiotherapy Equipments,
Occupational Therapy, Slimming Equipments**

Address : 315, RG Complex-1, Prashant Vihar,
Sector-14, Rohini, Delhi-110085

Phone No. : 011-47543164

Mobile No. : +91-9999303118, 9953888183

Email ID : info@accordmedicalproducts.com

Website : www.accordmedicalproducts.com

LASER THERAPY UNITS

AMP:03LA01: Computerized Cont. & Pulsed Lazer Therapy Unit (Super Deluxe)

FEATURES :

Versatile, user- friendly.

30 Inbuilt Memory Settings, enabling vast parameter storage.

Safety Lock to prevent accidental exposure.

2 Optional Hand Pieces 650 nm/ 850 nm (IR)

Highly user-friendly Parameter (Display) Panel, displaying Time, Frequency, Energy Level & Given Energy.

Compact build, economic usage.

TECHNICAL SPECIFICATIONS :

INPUT SUPPLY : 230V AC, 50 Hz

SIZE : L11.5" x W8" x H4.5"

WEIGHT : 2 Kg. Approx.

Programme : 00 Manual Programme

Number : 30 Empty Programmes for Memory Storage

Frequency : 25 Hz to 1500 Hz (It can be increased or decreased in multiple of 25 Hz and 100 Hz)

Energy : 3.0 m.J. (Min) To 5400 m.J. (Max)

Max. Power : 16 mW

Wave Length : 650 nm or 850 nm (I.R.)

Timer : 1 to 30 Minutes

STANDARD ACCESSORIES :

Laser protective goggles.

Hand piece (Probe) safety box.

Spare Key for safety lock.



AMP:03LA04: Lazer Therapy Unit(Probe Power :10Mw)(Economy)

TECHNICAL SPECIFICATIONS:

Mode: Continous and Pulse Mode

Max Probe Power: 10MW

Beam Surface: 0.025cm²

Laser Type: Semiconductor GaAs

Wave Leangth: 650nm

Power Density: 0.8w/cm²

Max Energy within 1 min 1 Joules



AMP:03LA03: Digital Laze

FEATURES :

TeraMed, STL-101 is very simple and easy to use : it has to be positioned so that the fix emission reaches the middle point of the part to be treated

TeraMed, STL-101 is completely automatic

TeraMed, STL-101 is successfully used in Dermatology, Aesthetic Medicine, Physiatrics, Rheumatology and pain therapy



TECHNICAL SPECIFICATIONS :

LASER SOURCE : Diode

MAIN POWER : 230 V AC, 50 Hz

DIMENSION : L130mm X W 36mm X H 45mm

WEIGHT : 28 Kg.

Scanner

Output power : 50mW

Wavelength : 650 nm

Automatic Scanner : Circle, Vertical, Atomic, Spot and combination

Radiation Time : 1~60 min

AMP:03LA05 Laser Therapy Unit(Probe Power:20mw)(Economy)



TECHNICAL SPECIFICATIONS:

Mode: Continuous and Pulse Mode

Max Probe Power: 20MW

Beam Surface: 0.025cm²

Laser Type: Semiconductor GaAs

Wave Leangth: 650nm

Power Density: 0.8w/cm²

Max Energy within 1 min 1 Joules

AMP:034902 Pathway™ MR-25 Dual Channel **EMG** Biofeedback System (with Vivid Tri-Color LED Display)



Continuous operation or work/rest prompts

Automatic data storage

Computer output

Pathway Utilities/Compliance software included for data download and real-time graphics

Hand-held portability

Carry case

Convenient 9 volt battery

Easy to operate

Audio and visual goal setting Include: Above Tone, Below Tone, Above Stim, Below Stim,

Maximum Display with Marker.

Includes: Pathway™ Utilities/Compliance PC Software and Computer Serial Interface Cable

ULTRA SOUND UNIT

AMP:03US01 Ultrasound Unit (MINI, Economy)



AMP:03US02 Ultrasound Unit(Digital) (Economy)



AMP:03US05 Ultrasound Unit(Heavy Duty LCD Model, Deluxe)



Technical specifications

Ultrasonic nominal Frequency	: 1 MHz, or 3 MHz
Ultrasonic Output	: 15Acoustic Maximum
Ultrasonic power density mode	
Operation	: 2.2 Watts/Cm2 Continuous and pulsed mode
Pulse ratios	: 1:1, 1:3, 1:6, 1:9
Display	: LCD Display with Programmable
Timer	: 0 to 99 minutes

AMP:03US04 Ultrasound

FEATURES :

Highly user- friendly Control Panel for Ultrasound Therapy in 24 different Conditions.

Pre- programmed Parameters (intensity, mode & time of treatment).

Manual setting option.

Compact design, Versatile features, Powerful performance.

TECHNICAL SPECIFICATIONS :

INPUT SUPPLY : 230V AC, 50 Hz

SIZE : L14" x W8" x H4.5"

WEIGHT : 3 Kg. Approx

Programme : 24 Preset Programmes

Frequency : 1 MHz + 10%

Output Mode : Continuous & Pulsed | Variable Pulsed

Frequency : 16Hz, 48Hz, 99Hz

Intensity : Continuous Mode 0 to 2 Watts/cm² | Pulsed Mode 0 to 3 Watts/cm²

Time : 1 to 60 Minutes

Buzzer : Audio buzzer at the end of treatment time

STANDARD ACCESSORIES :

•Sonic Gel



Interferential Therapy Unit

AMP:03IF01 Interferential Therapy Unit(Digital Pre. Prog + Russian Current,)(Su

FEATURES :

With Russian current.

Universal & Economic Unit.

Highly user friendly control panel.

Complete computerized unit.



TECHNICAL SPECIFICATIONS :

INPUT SUPPLY : 230 V AC, 50Hz

SIZE : L 11.5" x W 8" x H 4.5"

WEIGHT : 2.5 Kg. Approx

Selection

Russian Current : 2500Hz Carrier freq.

Medium Frequency : 4 KHz, 2 KHz

Base : 0 to 100 Hz

Sweep : 0 to 150 Hz

Form : Triangular 6/6 Sec | Trapezoidal 1/5 Sec |

Rectangular 1/1 Sec

Mode : (4pv) 45° Four pole vector Interferential | (4pv)

90° Four pole Vector Interferential | (4p) Four pole

Interferential Therapy | (2p) Two pole Bipolar

Interferential Therapy

Output Mode : Pads

AMP:03IF05 Interferential Therapy Unit(102Prg.LCD Model) (Deluxe)



Interferential therapy (ift) can be used for pain control. The Neural fibers are stimulated in the tissue where the maximum interference occurs and the two currents cross.

29Programes for IFT & 73 Programs fro MST & TENS.

Display : LCD

Frequency : 2&4 KHz

Sweep AMF : 0to 100Hz

Therapy Mode : 2P, 4P, 4PV

Balance : Available

Main Input : 200-230 VAC, 50Hz

Timer : 0 to 15 min. Digital

Base AMF : 6 to 12 Sec.

Spectrum Program : Liner, Trapezoidal & Non-Linear

Intensity : 0 to 100Ma (Adjustable)

Weight : 2.5Kg (Approx)

AMP:03IF02 Interfrential Therapy Unit(Digital 32Pre. Prog, Economy)



TECHNICAL SPECIFICATIONS :

Input Voltage: 220V/Ac50Hz

Display: FND Display

Intensity: 99mA

IFT Medium Frequency: 4KHz

IFT Muscle Exec.: 2kHz

Base Freq.: 0-150Hz

Sweep freq.: 0-100Hz

Therapy Mode: 4PIFT, 2PVector, 2P Modulation

Sweep Prg: Rec. Wave-1/1

Triangular Wave: 2/2

Trapezoidal Wave: 6/6

AMP:03IF03 Interfrential Therapy Unit (Digital 16Pre. Prog, Economy)



AMP:03IF06 Interfrential Therapy Unit (LCD, 32 Prg. Economy)



AMP:03IF08 Interfrential Therapy Unit (Portable)



Handy digital IFT and EMS device with the most effective preset programs designed for continuous pain relief, hassle free treatment which can be self administered.

Channel:- 2

Five Programs:- IFT with various frequency sweeps


Special Program:- Unique Sequencing of IFT & EMS

Interface:- Large LCD Display

Selectable :- 2 Pole Or 4 Pole applications

Power Source:- 9V Battery/ Adaptor

AMP:03IF07 Vaccum Therapy Unit for IFT(Economy)

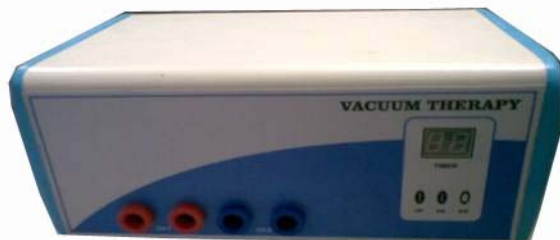


Short Wave Diathermy Units

AMP:03SW01 Short Wave Diathermy
(300 watt, Digital timer Economy)



AMP:03SW02 Short Wave Diathermy
(250 watt)Economy

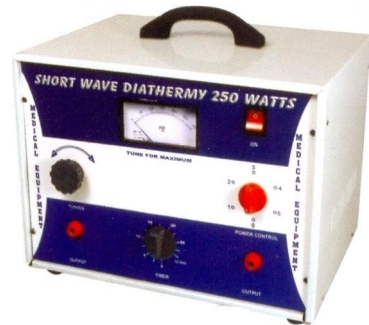

AMP:03SW03 Short Wave Diathermy
(500 watt, Digital timer, Pad Elec,) Optional: Sinus Elec, Disc Elec)Deluxe



AMP:03SW04 Short Wave Diathermy
(500 watt, Portable)Economy



AMP:03SW06 Short Wave Diathermy (500 watt) S.Deluxe
FEATURES :



5-STEP Intensity Control for smooth output.

Aesthetic design, with castor wheels (trolley not required)

30-minute Treatment Timer.

Automatic disconnection of HT Circuit on preset time lapse

TECHNICAL SPECIFICATIONS :

INPUT SUPPLY : 230 V AC, 50 Hz

SIZE : L20" X W13" X H 36"

WEIGHT : 40 Kg. Approx

OutPut Freq. : 27.12 MHz + 10%

OutPut : 500 Watts

STANDARD ACCESSORIES :

•One pair condenser pad & Cable

OPTIONAL ACCESSORIES :

•Disc electrodes

•Optional card update to PEME systems

TENS Units

AMP:03TN01 TENS
(2 Channel,
Portable Economy)



AMP:03TN02 TENS
(2
Channel, Imported)



AMP:03TN04 TENS
(4 Channel,
Digital Timer, Economy)




AMP:03TN02A Portable Analog TENS



Channels :- 2 with 4 Pads
Selection :- Pulse Width (40-260 μ S), Pulse Rate (2-150 Hz)
Modes :- Burst, Normal and Modulation
Power :- 9V Battery, 9V Adaptor
Features :- Low Battery indicator (LED) and Slide Cover
FDA USA 510(K) Cleared

AMP:03TN07 TENS
(6 Channel, Auto Mode,
Economy)

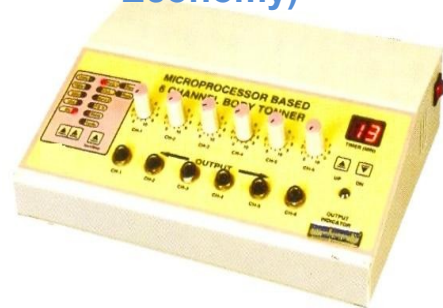

AMP:03TN03: TENS (4Channel, Deluxe)



Technical Specification

Frequency: 2Hz to 150 Hz
Pulse Width: 400 μ Sec
Patient Safety: Output is 100% Isolation
Therapy mode: All continuous, all burst, 2 continuous & 2 burst.
Output: 0-100V
Frequency indicator: Led indicator.
Digital Treatment Timer

AMP:03TN05: Computersied TENS (4Channel, with Digital Timer, Super Deluxe)

FEATURES :



Digital Timer with Buzzer

Rectangle- shaped Biphasic Pulse

Mode : Five types of Mode (Continuous, Fast Burst, Slow Burst, Ramp Burst and Sweep)

TECHNICAL SPECIFICATION :

·**Input Supply :** 230 V AC, 50 Hz

·**Size :** L 11.5" X W 8" X H 4.5"

·**Weight :** 1.5 Kg. (Approx)

·**Output Voltage :** 0 to 50 V with load of 500 ohm / 0 to 110 V without load

·**Pulse Width :** 50 to 350 Micro sec.

·**Frequency :** 16 to 200 Hz

STANDARD ACCESSORIES :

·Four Pairs Rubber electrodes , Four pairs Velcro Straps, One Gel Bottle, One Mains Cord

Muscle Stimulator

AMP:03MS01: Muscle Stimulator (Portable, Economic)



AMP:03MS06: Muscle Stimulator (Diagnostic, Economy)



AMP:03MS02: MINI Muscle Stimulator (Deluxe)



FEATURES :

Easy portable treatment

Three output selections (Tens, Interrupted Galvanic, Surge Faradic)

Good quality cover.

INPUT SUPPLY : 230 V AC, 50Hz

SIZE : L 6" x W 5.5" x H 2"

WEIGHT : 1.5 Kg. Approx

TNS

·Pulse width : 400 micro sec.

·Frequency : 50Hz

IG Pulse width : 100 mili sec.

·Frequency : 1 per sec.

SF : 1 to 9 sec. duration

STANDARD ACCESSORIES : One pair plate electrode, One pair Velcro Straps. One point electrode

AMP:03MS03: Muscle Stimulator (Di

FEATURES :

Highly user-friendly Control Panel

Universal & Economical Unit

Easy-to-operate design

SDC Facility available

Five types of Currents (Galvanic, Interrupted Galvanic, Faradic, Surge Faradic & TENS)



TECHNICAL SPECIFICATIONS :

INPUT SUPPLY : 230 V AC, 50Hz

SIZE : L 11.5" x W 8" x H 4.5"

WEIGHT : 2.3 Kg. Approx

Input Current : 58 mA

Output Current : 60 mA

Galvanic Output : Low ripple D.C. continuous adjustable from 0 to 145 volts

Interrupted Galvanic

Frequency : 1, 2, 3, 5, 10, 30, 50, 60, 100 Hz

Pulse Width : 0.01, 0.03, .10, .30, .40, 1, 3, 10, 30, 100, 300 Milli Seconds

Faradic : 1 Milli Second pulsed width and 60 Hz Frequency. Amplitude continuously adjustable from 0 to 145 volts

Surge Faradic

Contraction Period : Contraction period adjustable from 1 to 9 seconds

Rest Period : Rest period adjustable from 1 to 9 seconds

TENS : Pulse width 400 micro sec. & frequency 50 Hz

STANDARD ACCESSORIES : One pair plate electrode, One pair Velcro Straps, One point electrode

AMP:03MS04: Muscle Stimulator (Diagnostic, Deluxe)



Technical specifications

Plain Galvanic: It is for the stimulation of denervated muscle ion transfer.

Maximum output 75 MA.

Interrupted Galvanic Current: This is used for diagnosis and the extent of muscle lesion.

It is used for diagnosis of muscle activities by plotting S-D Curve of the defective muscle.

S-D Curve can be plotted on a semi log graph sheet. Pulse width variations are in 10 steps from .01 to 300MS.

© Faradic Current : Pulses generated are of short duration (0.6ms) at the rate of 80 pulses/second.

Tetanic muscles contraction and muscle fatigue is achieved in short duration since the rate is 80 PPS

(D) Surged Faradic Current: Faradic Pulses are surged. Peak tatanizing is reached from 2 Seconds.

to 4 seconds. Rest period Is adjustable from 2to 5 seconds. Action corresponds similar to voluntary Contraction.

COMBO Therapy Units

AMP:03CM09: IFT/TENS/EMS Handy Stimulator



Channels:- 4 with 8 Pads
Modes :- TENS (Burst), EMS (Continuous, Russian) and IFT (Low, High, Wide)
Power :- 5.5 V Adaptor with rechargeable battery
Features :- Work on Battery upto 2 Hrs.
FDA USA 510(k) Cleared

AMP:03CM02: COMBO (US+MST+TENS Dual, Deluxe)



AMP:03CM03: LCD Display COMBO (IFT+MS+TENS Dual, Economy)



AMP:03CM04: COMBO (US+TENS Dual)(Economy)



AMP:03CM05: COMBO (MS+TENS Dual, Economy)



AMP:03CM06: COMBO (US+MS+TENS Dual, Economy)



AMP:03CM07: COMBO (Digital IFT+MS+TENS Dual)(Economy)



Chest Massagers

AMP:03MA01: Chest Massager (Heavy Duty, Imported)Deluxe




AMP:03MA02: Chest Massager (Imported, Economy)



SLIMMING EQUIPMENTS

Multi purpose Electrotherapy Appliance

AMP-03SL01 : Micro- Computerized Body shaper (Deluxe)

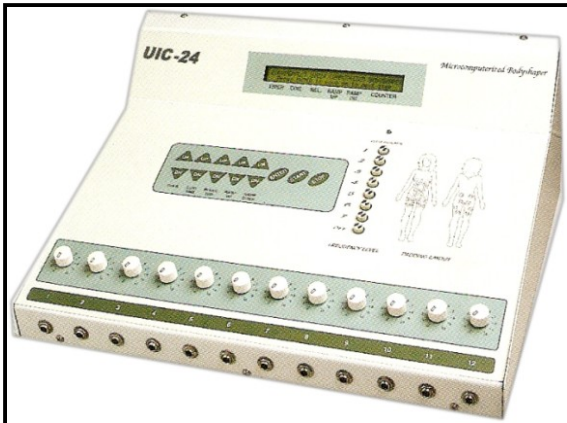

The most complete slimming machine for slimming centres. Easy to operate & variety of functions make this machine the most ideal machine for slimming centres. The instruments enables you to provide most versatile body treatments including obesity control, body contouring & firming treatments. All operations & parameters are monitored by easy to read LCD & controlled by computerized timing system.

Features

1. Adjustable contraction time of 1 to 30 seconds.
2. Adjustable Relaxation time of 1 to 30 seconds.
3. Adjustable ramp up time of 1 to 25 seconds.
4. Adjustable ramp down time to 1 to 25 seconds.
5. Adjustable timer for 1 to 60 seconds.
6. Counter to keep records of contractions.
7. Seven different frequency levels.
8. 24 pads to cover the entire body.
9. Enter key to input data of your choice.
10. Pause key to stop treatment in between.
11. LCD display.

AMP-03SL02 : Multipurpose Electrotherapy Appliance 24 Pads (Deluxe)



Used for slimming, Refirming, & anti-cellulite treatments. The most important body treatments are combined in this sleek appliance.

Features:-

1. Tapping, Blinking & stimulation -3 different frequency levels.
2. Continuous & Intermittent pulse mode selection.
3. 24 pads to cover the entire body –stomach. hips, thighs ,& arms.
4. Contraction & Relaxion time control knobs.
5. Digital timer to fix the time duration of treatment.

AMP-03SL03 : Multipurpose Electrotherapy appliance 20 Pads (Deluxe)



Easy to use , it has 10 currents outputs for treating the different parts of the body. All the important function of a professional slimming machines have been combined in this appliance for slimming, refirming & anti – cellulite treatments.

Features:

1. Adjustable contraction & Relaxion time settings in 3 pre – fixed time duration levels.
2. Tapping, stimulation ^ combination frequency option.
3. The trapping frequency can further be adjusted to 5 different speed levels.
4. 20 pads to cover the entire body –stomach, hips thighs & arms.
5. Digital timer to adjust the time duration of treatment.
6. Polarity change switch.
7. Start & pause key to start or pause the treatment in between.

AMP-03SL04 : Muscle Stimulation with comfort & Precision Deluxe




Designed primarily for body refirming this machine helps to restore the body's shapeliness and figure very effectively. It has 4 current outputs with 16 pads designed for toning the flabby muscles on the specify key areas like abdomen hips & thighs.

Features:

1. Specially designed for body toning.
2. 16 pads to cover the entire body.
3. Digital timer with display to fix the treatment time

AMP-03SL05 : Multipurpose electrotherapy appliance 36 Pads Delexe

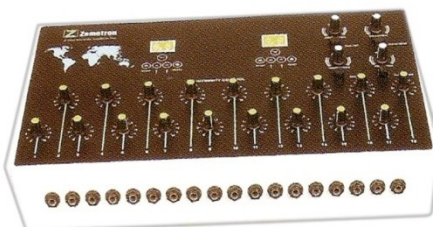


Used for slimming, refirming & anti –cellulite, treatments, this machine has been designed to treat two clients at one time. All the important controls are separated for the two clients being treated at one time. One client can be treated for refirming treatment whereas the other can be treated for anti-cellulite at the same time.

Features:


1. 18 channels with 36 pads .
2. 2 in 1 machine specially designed to treat two clients at one time .
3. It has separate controls for contraction & relaxation time.
4. Separate controls for frequency adjustment.
5. Separate timer with display for both clients.

Personal Slimming Machines

AMP-03SL08 : Slimming
Machines(16Pads)Deluxe

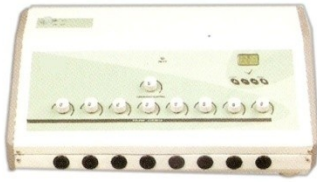

AMP-03SL06 : Slimming
Machines(16Pads) Deluxe

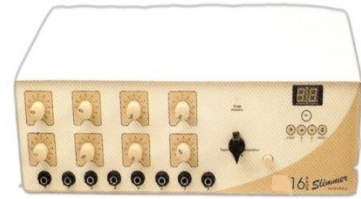

AMP-03SL13 : Slimming
Machines(12Ch.) Economy



AMP-03SL14 : Slimming
Machines(8Ch) Economy



Fat Analyser



AMP-03SL07 : Slimming
Machines(8Pads) eluxe



Deep Heat Therapy Units


AMP:03SL09 Deep Heat Therapy Units Deluxe



AMP:03SL17 Deep Heat Therapy Units Economy



AMP:03SL10 Deep Heat Therapy Units Deluxe

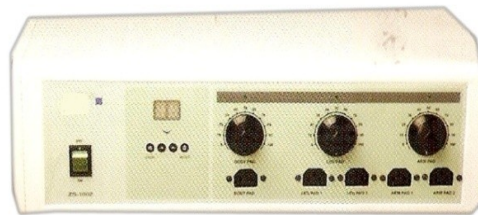


Vacuum Therapy Units

AMP:03SL11 Vacuum Therapy Units Deluxe



AMP:03SL20 Vacuum Therapy Units Economy

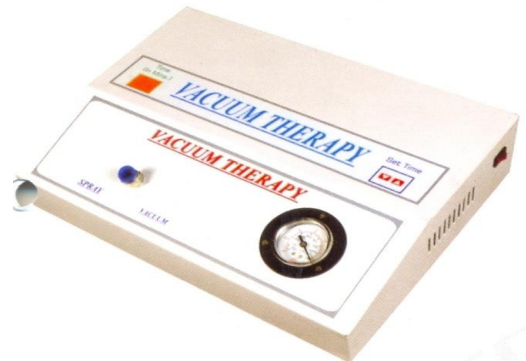

AMP-03SL12 : Professional Massage Apparatus Floor Model



With this powerful & durable massage apparatus you can give various high – demand services to your clients.

Accessories

1. Applicator for gentle massage of neck, shoulder, thighs & arms.
- 2. Applicator for rigorous massage of large body parts as hips buttocks & thighs.
- 3. Applicator for trigger point massage of body areas.
- 4. Applicator for deep massage of body areas like back, thighs, & buttocks.
- 5. Applicator of gentle massage of delicate & tender body parts.
- 6. Applicator for scalp massage.
- 7. Applicator for oil massage of neck, shoulder, & arms.